



# The EZ-Change® Attachment System

## Adaptation to Flexi-Overdenture® & EDS' AccessPost™ Overdenture



89 Leaning Street  
S. Hackensack, NJ 07606

1-800-22-FLEXI

CE APPROVED PRODUCTS  
Certificate No. CE 01420



Certificate No. FM 35005  
ISO 9001 / EN 46001




← **EZ-CHANGE® METAL KEEPER:**  
permanently cured into denture to hold nylon cap insert

← **EZ-CHANGE NYLON CAP INSERT:**  
easily removed, replaced when worn

This patented keeper and cap insert allows for "quick and easy" nylon cap replacement. When worn, the cap insert may be easily removed from the permanent metal keeper, requiring only seconds for replacement.

### OVERDENTURE SYSTEMS FEATURES & BENEFITS

## Flexi-Overdenture® & EDS' AccessPost™ Overdenture

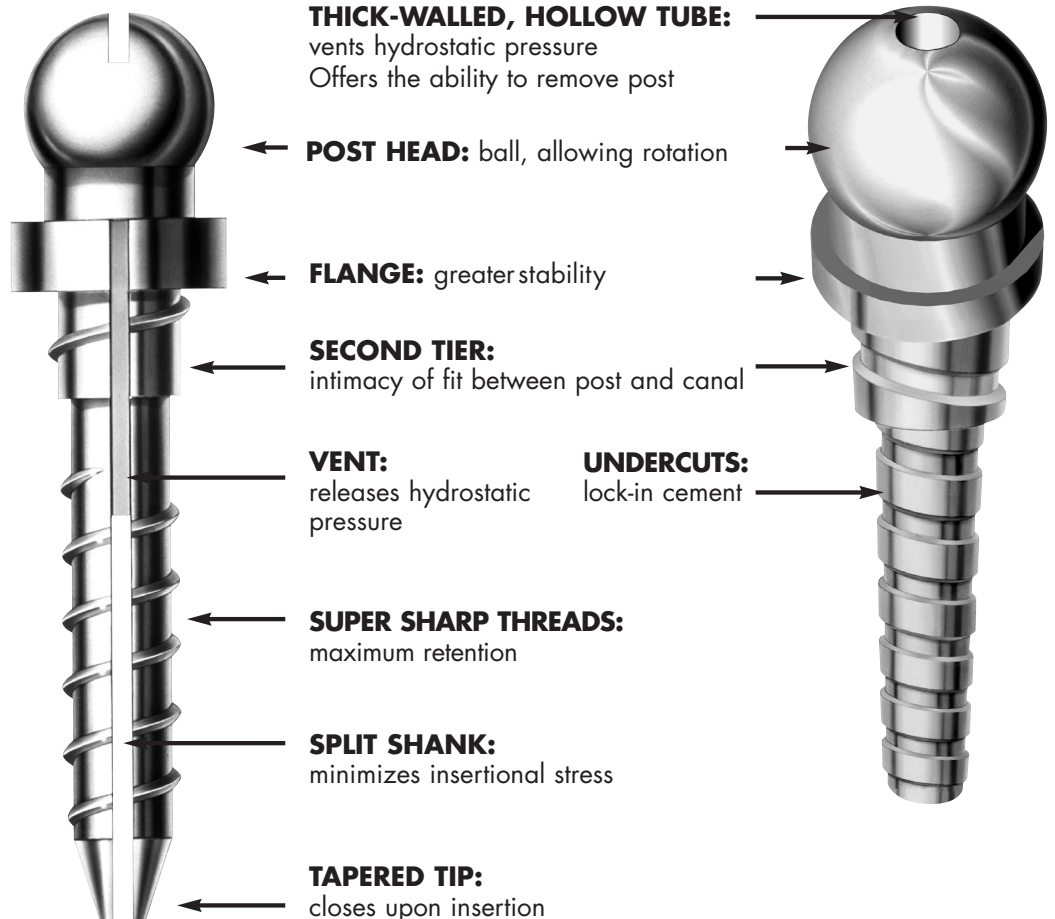

← **THICK-WALLED, HOLLOW TUBE:**  
vents hydrostatic pressure  
Offers the ability to remove post

← **POST HEAD:** ball, allowing rotation

← **FLANGE:** greater stability

← **SECOND TIER:**  
intimacy of fit between post and canal

← **VENT:**  
releases hydrostatic pressure


← **SUPER SHARP THREADS:**  
maximum retention

← **SPLIT SHANK:**  
minimizes insertional stress

← **TAPERED TIP:**  
closes upon insertion

← **UNDERCUTS:**  
lock-in cement


## EZ-CHANGE TECHNIQUE:

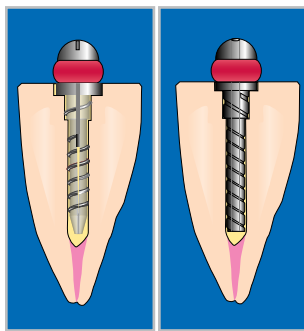


Flexi-Overdenture EDS' AccessPost Overdenture  
Overdenture post  
cemented in place.



Cap Insert is threaded  
into the metal keeper with  
2-pronged wrench.

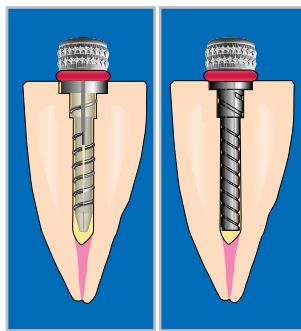

Flexi-Overdenture EDS' AccessPost Overdenture  
Cap Insert & Keeper  
placed on post.

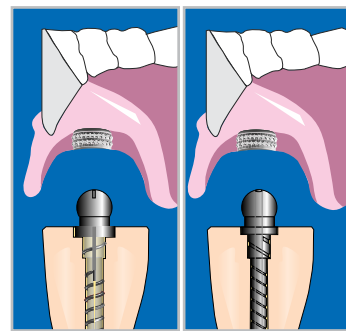

Flexi-Overdenture EDS' AccessPost Overdenture  
Cap Insert & Keeper  
incorporated within denture.

## EZ-CHANGE KITS AND THEIR CONTENTS

**Flexi-Overdenture EZ-Change Kits:**  
(2 keepers, 2 cap inserts, 2 rubber bands &  
EZ-Change® wrench)\* . . . . . **Cat. No. 250-00**  
(2 keepers, 2 cap inserts & 2 rubber bands)\*  
**Cat. No. 250-01**  
(6 nylon replacement cap inserts)  
**Cat. No. 250-02**  
(EZ-Change wrench) . . . . . **Cat. No. 250-03**

**EDS' AccessPost™ Overdenture EZ-Change Kits:**  
(2 keepers, 2 cap inserts, 2 rubber bands &  
EZ-Change® wrench)\* . . . . . **Cat. No. 550-00**  
(2 keepers, 2 cap inserts & 2 rubber bands)\*  
**Cat. No. 550-01**  
(6 nylon replacement cap inserts)  
**Cat. No. 550-02**  
(EZ-Change wrench) . . . . . **Cat. No. 550-03**


**\*COMPLETE INSTRUCTIONS INCLUDED**

## USING THE RIGHT PRODUCT IS ESSENTIAL


- FLEXI-POST®** - The patented split shank post is ideal where adequate coronal dentin exists. Maximum retention with minimal insertional stress.
- FLEXI-FLANGE®**-All the benefits of the patented split shank Flexi-Post while incorporating a flange to provide stability in situations where there is inadequate coronal dentin and excessive occlusal forces — anterior and posterior.
- FLEXI-OVERDENTURE®** - The patented split shank post with a ball and socket attachment that delivers the highest retention for full and partial dentures. Provides a simple, inexpensive overdenture at chairside.
- EZ-CHANGE®**- A patented keeper and cap system that allows for “quick & easy” nylon cap replacement. When worn, this cap may be easily removed, requiring only seconds for replacement. Available for use with **EDS' AccessPost™** Overdenture and Flexi-Overdenture systems.
- TI-CORE®** - Patented, reinforced, fluoridated core material. Cuts like dentin. No grooving or ditching. Available in original gray color or Ti-Core® Natural in Shade Vita® A3.
- FLEXI-FLOW Cem®** - A patented, reinforced, fluoridated composite cement with strength similar to dentin for superior post retention. Available in gray color or Flexi-Flow Natural® in Shade Vita® A3.
- EZ-FILL®** - The patented bi-directional spiral & epoxy root canal cement combined with a single point technique has been shown to create a seal equivalent to lateral condensation & thermoplastic gutta percha, both at a cost & time savings.
- EZ-FILL® SAFESIDER®**- A series of patented, non-circular, non-interrupted, flat-sided endodontic instruments. The patented EZ-Fill® SafeSider® system provides space for dentinal debris resulting in fast non-binding access to the apex.
- GLOVE'N CARE®** - A specially formulated, hypoallergenic hand cream for the care and relief of dry, irritated hands caused by wearing latex gloves. Glove'n Care is water based for the uncompromising maintenance of adhesive dentistry.




**INSTRUMENTATION  
& OBTURATION  
SYSTEMS**

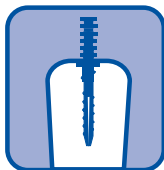

**POST  
SYSTEMS**

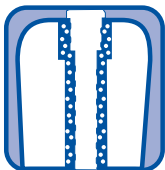

**COMPOSITE  
CEMENTS**

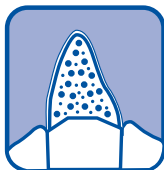

**CORE  
MATERIALS**



© 1/00 Essential Dental Systems, Inc.

89 Leuning Street, S. Hackensack, NJ 07606

**1-800-22-FLEXI**  
(1-800-223-5394)

TEL: (201) 487-9090 E-MAIL: eds@pipeline.com  
FAX: (201) 487-5120 http://www.edsdental.com

*All Essential Dental Systems Products Are  
Backed By A 30-Day Satisfaction Guarantee*