

Flexi-Flow Auto® E



Patented, lanthanide* reinforced, self curing, multi-purpose cements.

FEATURES

Titanium and lanthanide reinforced.

Auto mix, low viscosity.

Self curing formula.

Chemically bonds with composite core buildup material.

Radiopaque.

Insoluble.

Flexi-Flow Auto® E - Shade Vita® A2.

Compatible with all Ti-Core® reinforced core materials.

BENEFITS

Provides strength comparable to dentin.

Auto mixed and easily placed. Thoroughly coats dental tubules and post. Repeatable mix consistency.

Sets and hardens in 5 minutes, saving valuable chair time. Fully cures even at the bottom of the post-hole.

Saves you time. No removal of excess flash.

Clear x-ray viewing.

For lasting strength.

For aesthetic restorations.

Provides a homogeneous mass for solid core build-up.



Flexi-Flow Auto E - The Essential Link!

Eliminate what has always been the weakest link in any restorative system: the post-to-dentin interface. Since it surpasses the compressive strength of dentin, Flexi-Flow can be relied upon to maintain the strength found in natural tooth structure.

Now this rock solid reinforced flowable composite cement comes in a efficient delivery system that allows cement placement intra orally down the root canal!

**Flexi-Flow and Flexi-Flow Natural:
Available in original gel/gel consistency!**



* Data on file


Flexi-Flow

THE ESSENTIAL CEMENT.

Flexi-Flow's strength comes from varied particle and density size efficiently displacing oxygen from the matrices.

Flexi-Flow's high-density large particles and low-density small particles work together to form a highly compressive and efficient composite.

Flexi-Flow's matrix, reinforced by the cohesive strength of titanium, lanthanide, and silica keeps the post-to-dentin interface from breaking down.

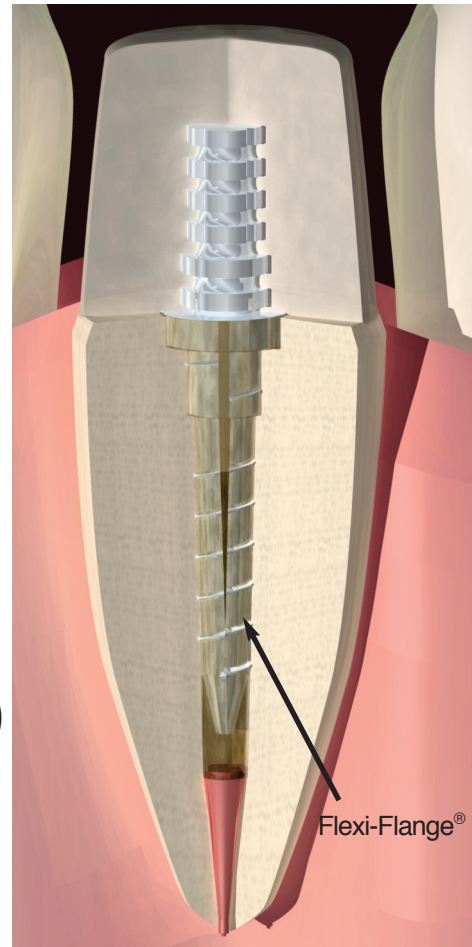

TECHNICAL FACTS:

Post Retention (w/ #1 Flexi-Post)	(985.2N (±117.9 lbs.) 221.5 lbs (±26.5 lbs.))
Diametral Tensile Strength	31.12 MPa (±4.27 MPa)
Compressive Strength	191 MPa (±22 MPa)
Radiopacity	Positive

ORDERING INFORMATION:

Flexi-Flow Auto E - Self Curing Reinforced Cement

Flexi-Flow Auto E Catalog Number 880-00

Includes:

2 auto mix syringes (each containing 4.5 gm base, 4.5 gm catalyst), 20 mixing tips, 20 intra oral tips.

Flexi-Flow & Flexi-Flow Natural

Flexi-Flow Catalog Number 850-00
Flexi-Flow Natural Catalog Number 860-00

Flexi-Flow Sets Include:

One 9g syringe of Flexi-Flow Base, One 9g syringe of Flexi-Flow Catalyst, 50 Disposable Spatulas, 50 Sheet Mixing Pad, and an Instruction Booklet

Flexi-Flow Auto Root Canal Intra-Oral Tips (ø 1.0mm OD)

20 Tips with Mixers (Fits all standard auto-mix double barrel syringes) Catalog Number 870-20

For our complete line of patented dental products please visit our website at: www.edsdental.com